




Aida®



SHORT STEM SYSTEM



Aida® - cementless short stem system

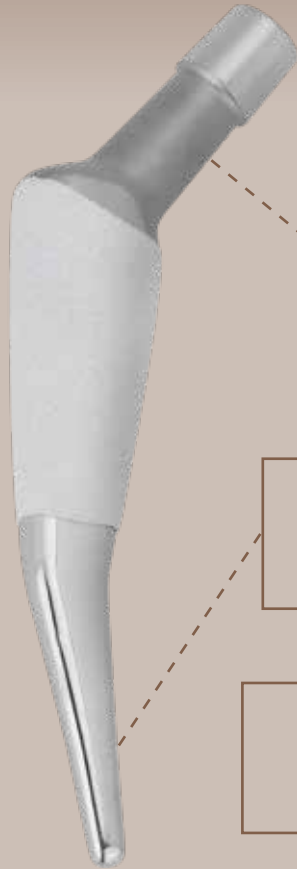

CCD angle:
130° (standard)
125° (lateralised)

material:
implatan® TiAl₆V₄

**uncoated distal tip for
intramedullary guidance**

proximal coating:
implaFix® Duo cpTi + HA

size spectrum:
nine sizes each standard and lateralised

coating thickness:
cpTi 250±50µm
HA 90±30µm



Aida

combination with ic-head	CoCrMo	Titanium with TiN- coating	BIOLOX® forte	BIOLOX® delta
Ø 22mm S, M, L	x			
Ø 28, 32, 36mm S, M, L, XL	x	x	x*	x**
Ø 28, 32, 36mm XXL, XXXL	x	x		
Ø 40mm S, M, L, XL				x

*BIOLOX® forte not available in XL

**BIOLOX® delta not available in Ø 28mm XL

metaphyseal fit

**stem design:
double-tapered shape,
trapezoidal cross-section**

partially neck preserving

**minimally invasive
instruments available**



implantcast GmbH
Lüneburger Schanze 26
21614 Buxtehude
Germany
phone: +49 4161 744-0
fax: +49 4161 744-200
e-mail: info@implantcast.de
internet: www.implantcast.de

CE 0482

